


Strip Punch Cutter

Product code [RL-SPC-A](#)

RYCOLAB

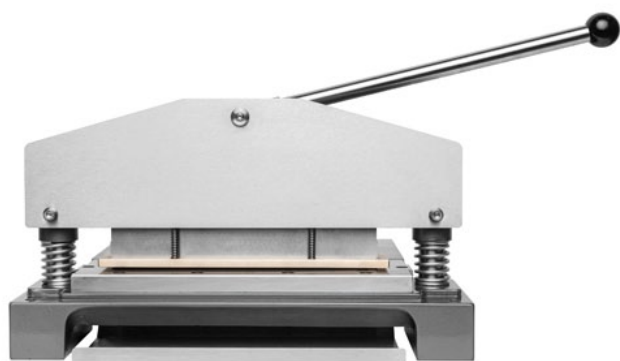

Usage

For cutting strips according to the »punch and die« principle.

Applicable standards

ISO 7263, DIN-EN ISO 7263, TAPPI T809, SCAN P27, PAPTAC D29, etc.

Description

Sturdy punch fully made from metal with feeder table and hand lever for easy operation. The punch and die cutting unit consists of an upper and lower blade and guarantees accurate sample strips. The sample material is placed on the feeder table and supported by the feeder edges placed under the punching unit. By pushing the hand lever down the sample is clamped automatically and a safety guard covers the punching unit. The samples strips fall through the lower part of the punching unit and are collected on a metal plate.

Dimensions

	Net	Gross
Width	600	800
Depth	500	600
Height [mm]	450	500
Weight [kg]	45	80